

HARMONIC ACE® Curved Shears

SPECIALTIES:

- General Surgery
- Colorectal
- Gynecology
- Urology
- Bariatric



DEVICES							
CODE	DESCRIPTION	GRIP	BLADE	SHAFT DIAMETER	LENGTH	COMPATIBLE WITH	QTY
ACE23E	HARMONIC ACE® 23	Ergonomic	Curved	5mm	23cm	HP054	6
ACE36E	HARMONIC ACE® 36	Ergonomic	Curved	5mm	36cm	HP054	6
ACE45E	HARMONIC ACE® 45*	Ergonomic	Curved	5mm	45cm	HP054	6
ACE14S	HARMONIC ACE® 14S	Scissors	Curved	5mm	14cm	HP054	6
CAPITAL							
CODE	DESCRIPTION	QTY	COMPATIBLE WITH				
GEN04	Generator 300 (110V)	1	ACE45E, ACE36E, ACE23E, ACE14S				
CRT01	Generator Cart	1	N/A				
FSW01	Footswitch and Cable	1	ACE45E, ACE36E, ACE23E, ACE14S				
HP054	HARMONIC® Hand Piece	1	ACE45E, ACE36E, ACE23E, ACE14S				
HST01	Sterilization Tray	1	N/A				
HST02	Test Tip	1	ACE45E, ACE36E, ACE23E, ACE14S				
TLB01	Torque Lock Blade Wrench	2	ACE45E, ACE36E, ACE23E, ACE14S				

*Coming soon

For more information, or to arrange a product evaluation, call your Ethicon Endo-Surgery Representative at 1-800-USE-ENDO.

Precision, multifunctionality and minimal thermal tissue damage. **That's positive energy.**



HARMONIC ACE® Curved Shears

Precision, multifunctionality and minimal thermal tissue damage.

That's positive energy.

ORDERING




HARMONIC ACE® Curved Shears

For laparoscopic and open procedures

Design refined by experience

Input from hundreds of surgeons along with continual refinements applied over the years have combined to produce a surgical device with maximum versatility and precision.

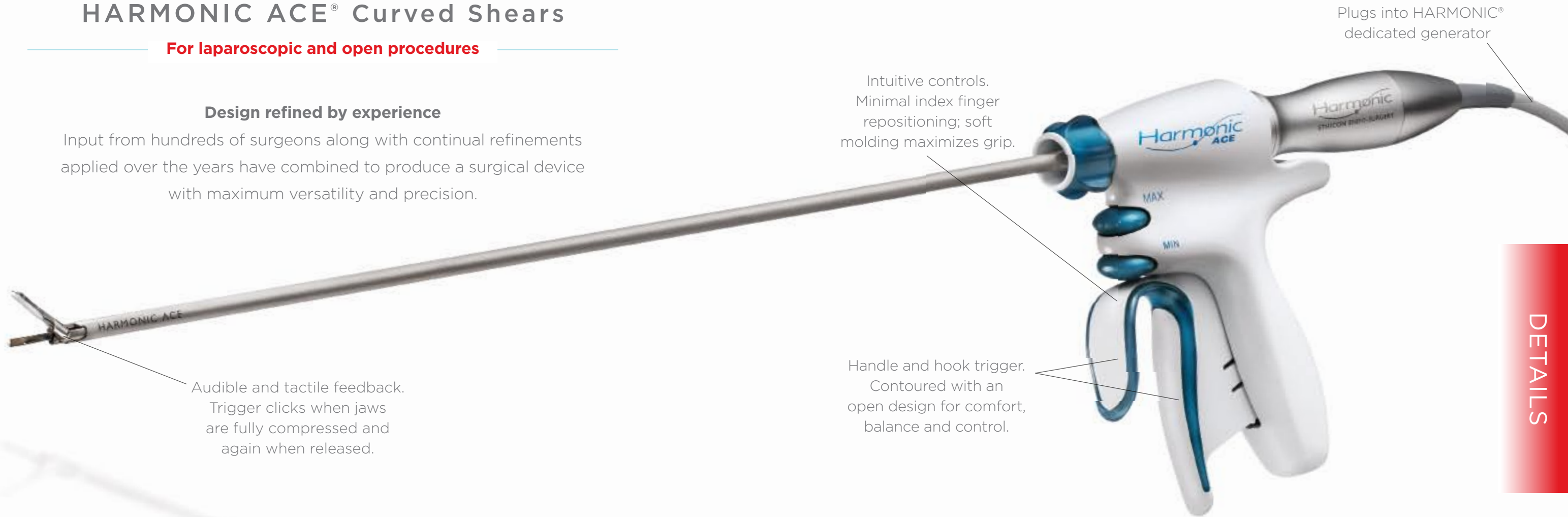


Plugs into HARMONIC® dedicated generator

Intuitive controls. Minimal index finger repositioning; soft molding maximizes grip.

Handle and hook trigger. Contoured with an open design for comfort, balance and control.

Audible and tactile feedback. Trigger clicks when jaws are fully compressed and again when released.

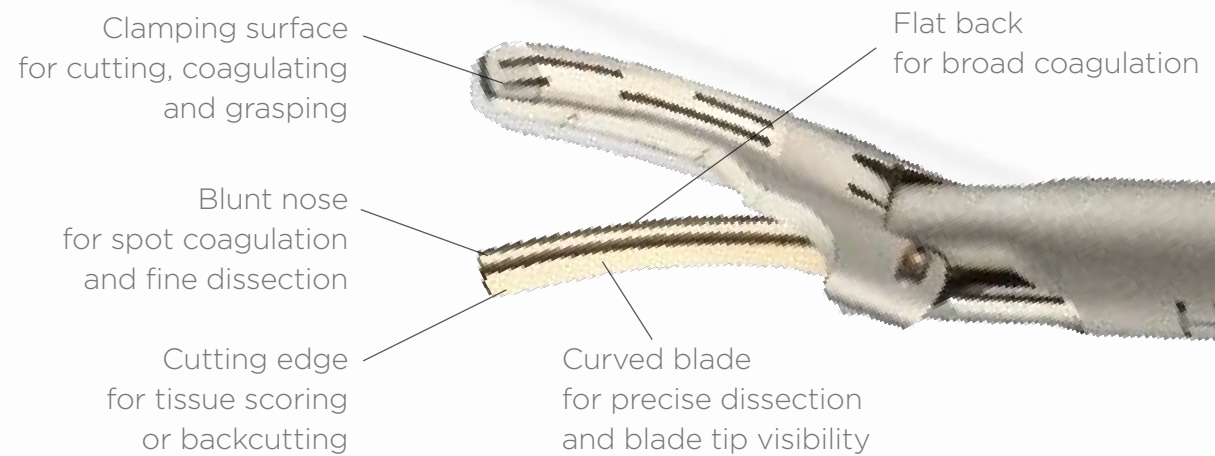

Clamping surface for cutting, coagulating and grasping

Flat back for broad coagulation

Blunt nose for spot coagulation and fine dissection

Cutting edge for tissue scoring or backcutting

Curved blade for precise dissection and blade tip visibility

DETAILS

PRECISE	Minimal lateral thermal tissue damage for safer dissection near vital structures
EFFICIENT	Moves through tissue quickly while maintaining hemostasis
VERSATILE	Precisely cuts, coagulates, grasps and dissects - without instrument exchanges
RELIABLE	Allows you to seal vessels ≤ 5mm, as well as lymphatics
ERGONOMIC	Comfortable grip and intuitive controls